


TESTING INSTRUMENTS FOR
A MEASURABLE DIFFERENCE...

IDM
instruments

CARTON CREASE PROOFER

MODEL: IDM-C0053-M1

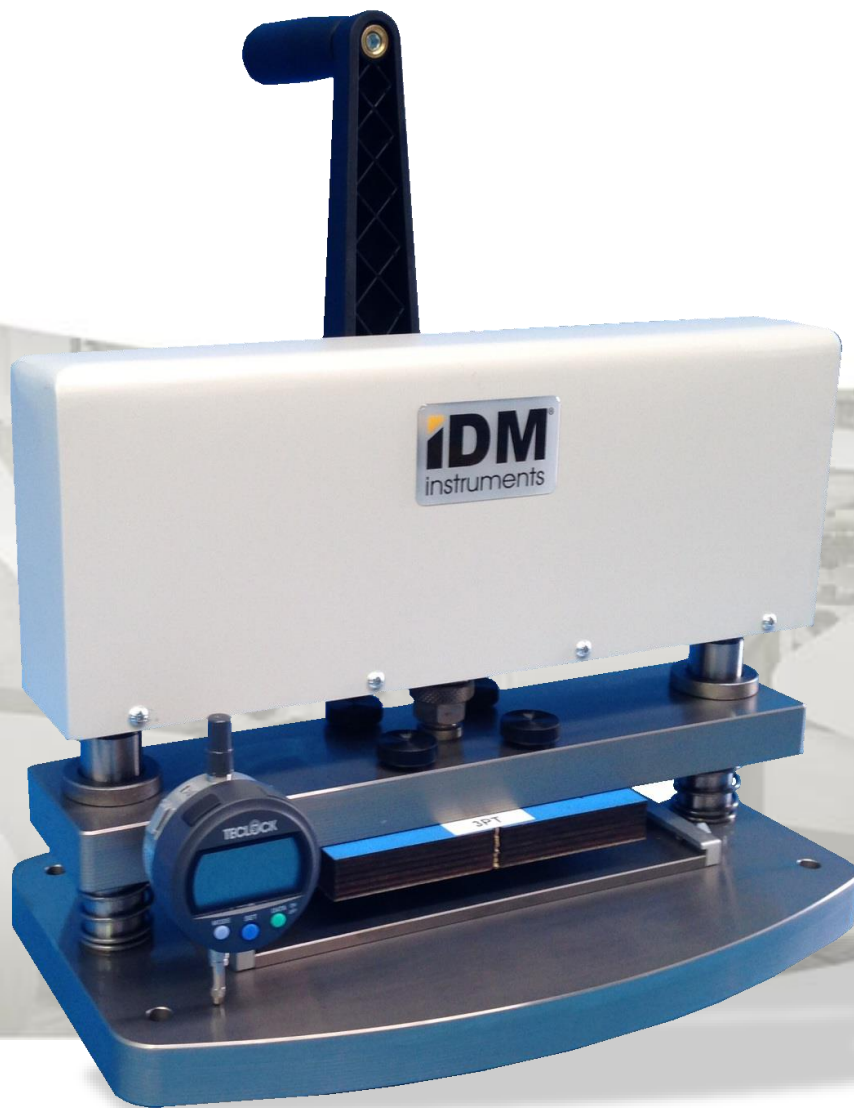


INTRODUCTION

This user-friendly instrument is commonly utilised in ink laboratories to assist the development of flexible coatings that resist cracking on commercial cartons. It allows the operator to quickly produce quality carton board crease samples without committing machine time.


The unit comes complete with robust rules & dies which are of the same quality of those used with a full scale cut and crease machine.

Result Guide Pictures


Beside main raised rib

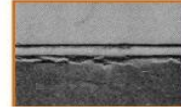

Across grain crumpled rib

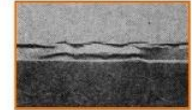

Along grain crumpled rib

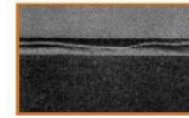

FEATURES

- Substrate: Suitable for any substrate 300 - 1000µm
- Creasing Rules: 2 PT, 3 PT, 4 PT
- 18 Crease Channels: 1 - 4mm
- Crease Depth: Variable

APPLICATIONS

- Carton Board
- Ink / Printing

BENEFITS

- Control Substrate Quality
- Determine Ink & Coating
- Flexibility

DIMENSIONS

Instrument:

- H: 350 mm
- W: 390 mm
- D: 200 mm
- Weight: 25 kg

Packaged:

- H: 500 mm
- W: 500 mm
- D: 600 mm
- Weight: 45 kg

WARRANTY & CALIBRATION

- 1 year Warranty
- Our **Preventive Maintenance and Calibration (PM&C) program** has been designed to make the maintenance and calibration of your valuable testing equipment more cost effective by preventing breakdowns and downtime by regular calibration, service and replacement of defective parts. Talk to us about this today.

RELATED ITEMS

- Crease & Stiffness Tester – IDM-C0039
- Crease & Stiffness Cutter – IDM-C0016

